


Klarity[®] UNO[™]



Precision
in Every
Position

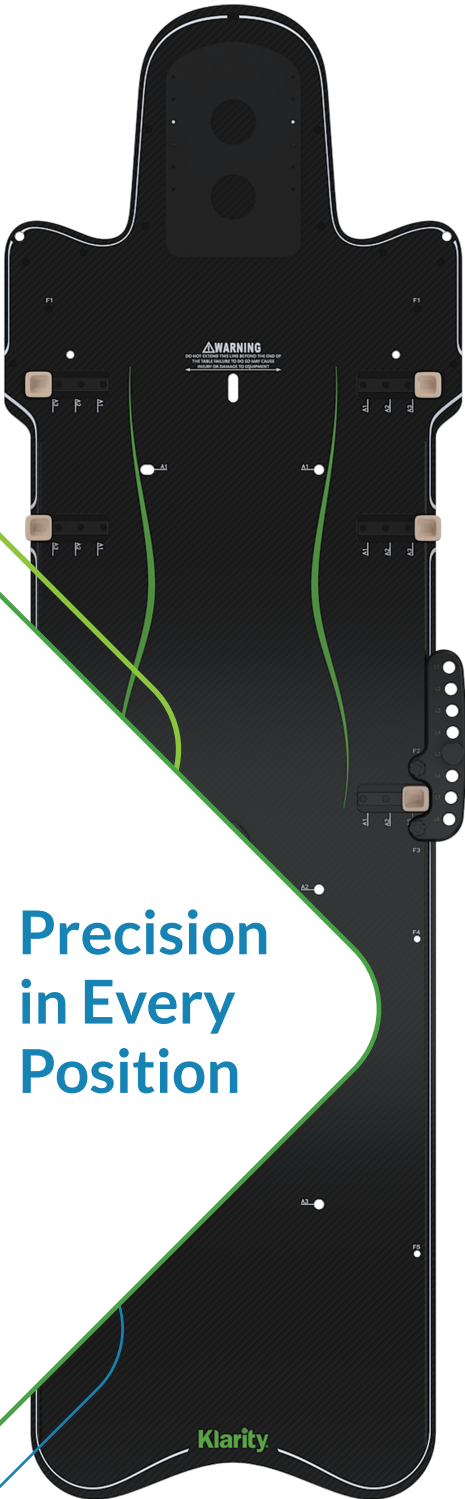

KLARITY UNO™

The Klarity Uno™ positioning system is highly adjustable and customizable and is suitable for all treatment types.

The Uno™ Baseplate comes as a set with a headrest adaptor, right and left handgrips, a groin lock, and can be equipped with standard S-type masks as well as Klarity Uno™ Breast, Torso and Pelvic Masks.

Additional configurations are possible, including a Tilting Board, Wingboard, Overhead Arm Positioner, SRS Adaptor, SBRT Set, Knee Positioner, Foot Positioner, and Prone Breast and Belly Systems. A storage cart and a tool accessory box are also available.


The Uno™ System is compatible with additional positioning and immobilization tools, including thermoplastic masks, Klarity Cushions™ and other head supports, vacuum bags, and more. Please see our website for a full list of compatible accessories at www.klaritymedical.com/uno.




KLARITY UNO™

The Uno™ is available in three sizes.

Standard

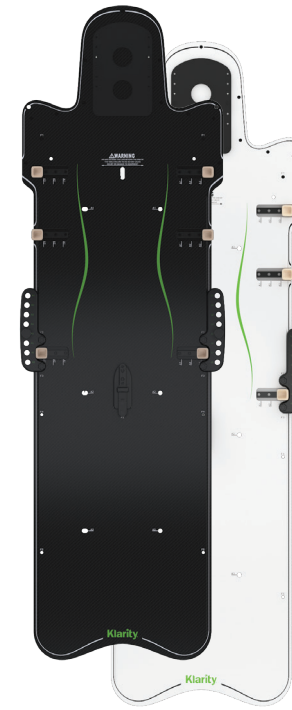

Narrow

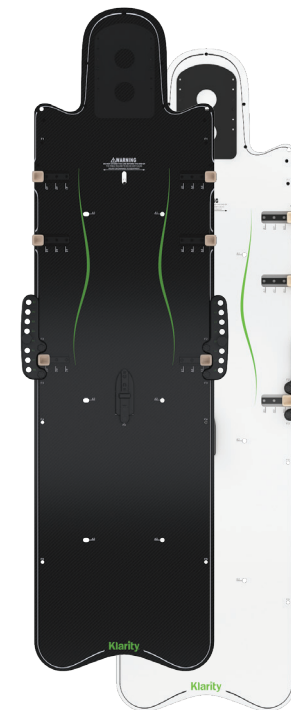

HexaPod



R616-SCFA	R616-SKFA
Carbon Fiber	Kevlar
199.5cm (length) x 61cm (shoulder width) x 54cm (body width)	199.5cm (length) x 61cm (shoulder width) x 54cm (body width)
Max Load 250kg	Max Load 250kg
5.6kg	

R616-NCF	R616-NKF
Carbon Fiber	Kevlar
199.5cm (length) x 54cm (shoulder and body width)	199.5cm (length) x 54cm (shoulder and body width)
Max Load 250kg	Max Load 250kg

R616-SCF1A	R616-SKF1A
Carbon Fiber	Kevlar
199.5cm (length) x 61cm (shoulder width) x 47.8cm (body width)	199.5cm (length) x 61cm (shoulder width) x 47.8cm (body width)
Max Load 250kg	Max Load 250kg

COMPATIBLE

All S-type Head Masks
All S-type Head & Shoulder Masks
Klarity S-type Pediatric Masks
Uno™ Breast, Torso, and Pelvic Masks
Silverman Head Support

COMPATIBLE








All S-type Head Masks
Uno™ S-type Head & Shoulder Masks (RLG460-2436S5T, R460-2436S5T)
Uno™ Breast, Torso, and Pelvic Masks
Silverman Head Support

COMPATIBLE

All S-type Head Masks
All S-type Head & Shoulder Masks
Klarity S-type Pediatric Masks
Uno™ Breast, Torso, and Pelvic Masks
Silverman Head Support

UNO™ ALL IN ONE SYSTEM

AVAILABLE COMPONENTS	PRODUCT ID	INCLUDED	UNO™ COMPATIBILITY
	R616-PC2 Wingboard (Carbon Fiber) R616-PK2 MR Wingboard (Kevlar) Compatible with the Athena™ Supine Breast System Can also be placed directly on the couch top	Wingboard T-bar Handles U-grip Handles Mask Adaptor Head Support	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD
	R616-PC1 Overhead Arm Positioner (Carbon Fiber)	Overhead Arm Positioner Baseplate 2 Arm Cups 2 Wrist Cups 2 Hand Grips	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD
	R616-PC3 Tilting Board (Carbon Fiber) R616-PK3 MR Tilting Board (Kevlar) Tilting Angles: 5°, 10°, 15°, 20°	Tilting Board Base Cushion 2 Lower Handgrips 2 Arm Cups (R616-PC3 only) 2 Wrist Cups 2 Upper Hand Grips Upper Indexing Bar Lower Indexing Bar	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD Requires one supplemental R616-PDCF-S4 Indexing Bar
	R616-PDCF-S3 Uno™ Tilting Board Anterior Indexing Bar Indexes supplemental head or head and shoulder overlays on top of the Tilting Board	Tilting Board Anterior Indexing Bar	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD
	R616-L1 Knee Positioner Includes 2 elevation blocks	Knee Positioner Baseplate Knee Riser 2 Elevation Blocks	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD
	R616-L3 Foot Positioner	Foot Positioner Baseplate Foot Holder	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD

AVAILABLE COMPONENTS	PRODUCT ID	INCLUDED	UNO™ COMPATIBILITY
	R616-AIR Compression Belt	Compression Belt Hand Pump	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD
	R616-PC4 Belly Bridge	Belly Bridge Belly Paddle Inflatable Belly Paddle Hand Pump	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD
	R616-SCF2-2 Knee Bridge	Knee Bridge	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD
	R7131-5NLTC Knee Bridge Vacuum Bag 28cm x 80cm 5-liter fill	Knee Bridge Vacuum Bag	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD
	R7644-32.5BC SBRT Vacuum Bag 136.5cm x 90cm 32.5-liter fill	Vacuum Bag with Two Chambers and Indexing Bar	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD Requires one supplemental RE-616-CF1 Indexing Bar
	RE-616CF Uno™ Lower Indexing Bar Attaches to the top of the Uno™ for Vacuum Bag Indexing	Lower Indexing Bar	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD
	RE-616CF1 Uno™ Hexapod Lower Indexing Bar Attaches to the top of the Uno™ Hexapod Baseplate for Vacuum Bag Indexing	Lower Hexapod Indexing Bar	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD

AVAILABLE COMPONENTS	PRODUCT ID	INCLUDED	UNO™ COMPATIBILITY
	S-20CF SRS Adaptor (Carbon Fiber) S-20A MR SRS Adaptor (Plastic)	SRS Adaptor	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD
	R616-BE5 Prone Breast Cushion Set Can also be placed directly on the couch top	Upper Prone Cushion Tilting Prone Head Support Right and Left Upper Hand Grips Flippable Breast Cushion Lower Prone Cushion with Removable Belly Hole Cushion Leg Bolster Upper Prone Cushion Indexing Bar Lower Prone Cushion Indexing Bar	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD <small>Requires two supplemental R616-AE5-SP1 Indexing Bars</small>
	R616-AE5 Prone Belly Cushion Set Can also be placed directly on the couch top	Upper Prone Cushion with Removable Port Cushions Tilting Prone Head Support Right and Left Upper Hand Grips Lower Prone Cushion with Groin Lock and Removable Pelvic Hole Cushion Leg Bolster Upper Prone Cushion Indexing Bar Lower Prone Cushion Indexing Bar	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD <small>Requires two supplemental R616-AE5-SP1 Indexing Bars</small>
	R616-BT1 Uno™ SBRT Set *Available in select markets only	Wingboard Belly Bridge Belly Paddle Inflatable Belly Paddle Hand Pump Knee Bridge Knee Bridge Vacuum Bag Knee Positioner Foot Positioner SBRT Vacuum Bag	<input checked="" type="checkbox"/> STANDARD <input checked="" type="checkbox"/> NARROW <input checked="" type="checkbox"/> HEXAPOD
	R616-CAB Uno™ Storage Cart	Uno™ Storage Cart	N/A

SBRT

The SBRT package offers unmatched versatility with arches, knee positioner, and foot positioner that can all be indexed anywhere along the baseplate. Three abdominal compression variations offer enhanced conformability.



SUPINE BREAST

The Uno™ features comprehensive Supine Breast positioning options, with a Wingboard, Overhead Arm Positioner, and a Tilting Board with elevation options that can be equipped with arm cups or Wingboard. Breast masks are also available. All supine breast set-ups can be used with lower hand grips for asymmetrical arm positioning.



PRONE BREAST AND PRONE BELLY

The Prone Breast and Prone Belly Cushion Sets can be used independently on a couch top. The Uno™ Prone Breast Cushions can be easily flipped over for Right/Left Positioning.




HEAD, NECK, AND TORSO



The Uno™ offers streamlined positioning for all Head, Neck, and Torso treatments. Masks click directly into the baseplate and the Uno™ is compatible with all standard positioning tools and accessories.


Compatible with S-type Masks and Uno™ Torso Masks

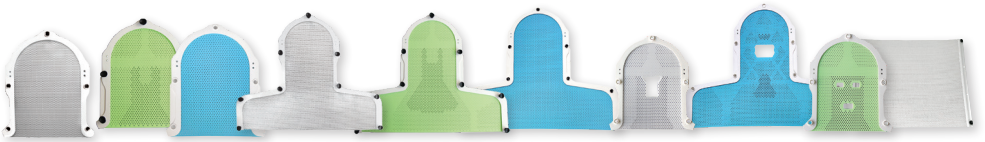

SRS

The SRS Adaptor connects easily and can be used with Klarity SRS Mask Sets. The two-piece mask systems allow for both anterior and posterior support while providing a 0-4mm range of adjustment.

Compatible with Klarity SRS Masks

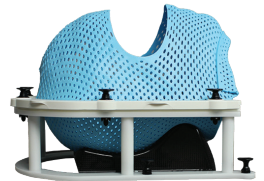

[KLARITYMEDICAL.COM](https://www.klaritymedical.com)